



BIOMÉRIEUX

VIDAS[®] B·R·A·H·M·S PCT[™]

Helping you manage
antibiotic therapy in LRTI*



PIONEERING DIAGNOSTICS

Did you know?



LRTI* account for **10%** of the worldwide burden of morbidity and mortality¹



75% of all antibiotic doses worldwide are prescribed for acute respiratory tract infections in spite of mainly viral cause¹



50% of antibiotics prescribed in the US for acute respiratory conditions are **unnecessary**²

VIDAS[®] B·R·A·H·M·S PCT[™]

Helping you manage antibiotic therapy in LRTI*

BECAUSE IT MAKES SENSE ON VIDAS[®]

- **On-demand** testing adapted to Primary Care and Emergency settings
- Rapid result: **only 20 minutes**
- Automated test on **easy-to-use** benchtop instrument
- **Cost-efficient** solution: 1 patient = 1 test


PRESCRIBE ANTIBIOTICS – OR NOT – WITH CONFIDENCE

To reduce antibiotic exposure and side effects and improve survival, you can count on VIDAS[®] B·R·A·H·M·S PCT[™]. It has been **validated in a clinical setting** for the management of antibiotic therapy in LRTI* patients¹².

Numerous international, multi-center, randomized control studies have shown that **PCT-guided antibiotic therapy in LRTI* is safe and effective**^{13,14}.

Significant reduction was observed in the PCT groups vs. the control groups in:

- **antibiotic initiation:** 70% vs 86%
- **antibiotic treatment duration:** 8 days vs 9.4 days
- **total antibiotic exposure:** 5.7 days vs 8.1 days
- **antibiotic-related side effects:** 16% vs 22%
- **mortality at 30 days:** 9% vs 10%

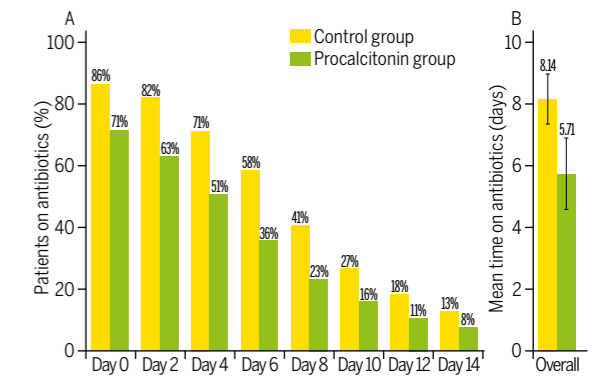

PCT reduces antibiotic exposure and treatment duration¹⁴
 A - Proportion of patients on antibiotics
 B - Mean duration of antibiotic use

The Challenge of LRTI* Management: Does Your Patient Really Need Antibiotics?

Many different forms:

- Acute bronchitis
- Community-Acquired Pneumonia (CAP)
- Acute exacerbation of chronic obstructive pulmonary disease (COPD)

Non-specific clinical characteristics:

- Cough
- Sputum
- Fever
- Shortness of breath

In the past, this challenge was met with empiric antibiotic therapy, resulting in:

- Antibiotic misuse
- Inappropriate initiation and prolonged use
- Antibiotic-associated infections such as *Clostridium difficile*, and other adverse effects^{3,4}
- Increasing antimicrobial resistance

Today, a valuable diagnostic tool can help you decide if antibiotics are warranted and optimize treatment duration.


How does Procalcitonin (PCT) Provide Critical Information?

A host biomarker specific for severe bacterial infection⁷⁻¹⁰

PCT is produced by numerous organs at a cellular level after bacterial pro-inflammatory stimulation^{5,6}. Its expression is inhibited by viral infections.

- Rises within 3-6 hours from bacterial insult
- Daily decrease of 50% as infection is resolved
- Proven to aid in monitoring patient response to antibiotic therapy¹¹

Monitoring antibiotic treatment with serial PCT measurements
 Adapted from Brunkhorst FM.⁹

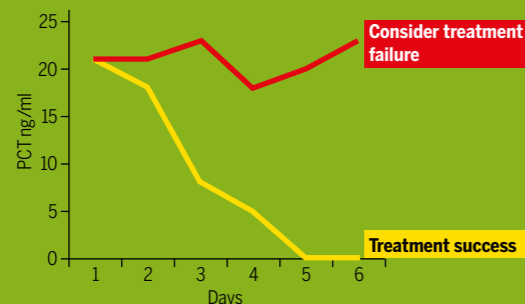

* LRTI: Lower Respiratory Tract Infections.



Give antibiotics to the right patients, at the right time, for the right duration

VIDAS® B·R·A·H·M·S PCT™ HELPS YOU IMPROVE ANTIBIOTIC STEWARDSHIP

START ANTIBIOTICS OR NOT?

Recommendations to start antibiotics¹²:

PCT value	< 0.1 ng/mL	0.1-0.25 ng/mL	0.26-0.5 ng/mL	> 0.5 ng/mL
Bacterial infection**	Absent	Unlikely	Possible	Suggestive
Initiation of antibiotic therapy	Strongly discouraged	Discouraged	Recommended	Strongly recommended

WHEN TO STOP ANTIBIOTICS?

Testing PCT at regular intervals*** combined with clinical assessments can support decision-making on antibiotic discontinuation for patients with LRTI*:

If PCT ≤ 0.25 ng/mL ► Discontinue antibiotics¹²

For cases of sepsis in the ICU:


If PCT < 0.5 ng/mL (or PCT drop ≥ 80%) ► Discontinue antibiotics¹¹

Consider continuing antibiotics in presence of clinical instability or disease progression.

If PCT remains high, consider treatment failure.

SUPPORTING HOSPITAL COST OPTIMIZATION

Several health economics studies have shown that PCT-guided antibiotic strategies in LRTI* are **cost-effective in different settings**: primary care, intensive care units and emergency departments.

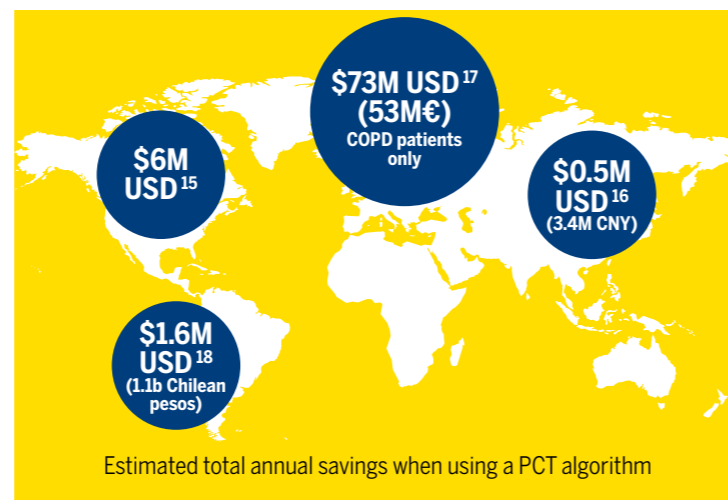

“Use of sensitive procalcitonin measurements in clinical algorithms can reduce antimicrobial overuse, decreasing the risk of side effects and controlling emerging bacterial multiresistance.”

Schuetz P, et al. Overview of procalcitonin assays and procalcitonin-guided protocols for the management of patients with infections and sepsis. Expert Review of Molecular Diagnostics 2017;17(6):593-601.



“PCT may be safely incorporated into treatment algorithms for children with CAP to reduce antibiotic administration and duration.”

Stockman et al. Procalcitonin Accurately Identifies Hospitalized Children with Low Risk of Bacterial Community-Acquired Pneumonia. Journal of The Pediatric Infectious Diseases Society 2017 Feb 3.doi:10.1093.



“Given the prevalence of COPD and the duration of illness, a reduction in antibiotic prescriptions for the treatment of exacerbations could have a tremendous impact on the overall economic burden of the disease under current budget constraints. In addition, the controlled prescription of antibiotics decreases selective pressure for the emergence of bacterial resistance.”



Stolz D, et al. Antibiotic treatment of exacerbations of COPD: a randomized, controlled trial comparing procalcitonin-guidance with standard therapy. Chest 2007;131(1):9-19.



“Good compliance with the PCT algorithm is possible in real-life conditions but has to be reinforced to achieve optimal benefit.”

Albrich WC, et al. Effectiveness and Safety of Procalcitonin-Guided Antibiotic Therapy in Lower Respiratory Tract Infections in “Real Life”: An International, Multicenter Poststudy Survey (ProREAL). Arch Intern Med 2012;172(9):715-722.

Discover more about PCT and antimicrobial stewardship:

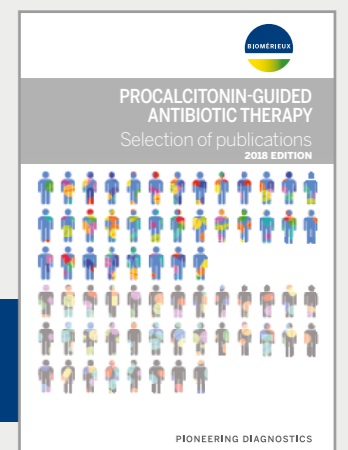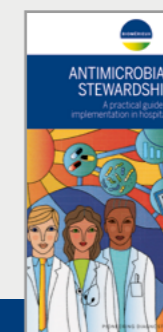

AVAILABLE ON

<http://www.biomerieux-diagnostics.com>
or from your local bioMérieux representative

** Adapted from VIDAS® B·R·A·H·M·S PCT™ package insert
*** Every other day or every 2 days, depending on the patient evolution.

AVAILABLE ON INSTRUMENTS OF THE VIDAS® FAMILY: VIDAS®, MINI VIDAS® AND VIDAS® 3



VIDAS® EMERGENCY PANEL

- VIDAS® B·R·A·H·M·S PCT™
- VIDAS® D-DIMER EXCLUSION™ II
- VIDAS® High sensitive Troponine I
- VIDAS® NT-proBNP2
- VIDAS® CK-MB
- VIDAS® Myoglobin
- VIDAS® Digoxin

VIDAS® B·R·A·H·M·S PCT™

Reference number	30450
Tests / kit	60
Time to result	20 minutes
Sample type	Serum, Plasma (Lithium heparinate)
Sample volume	200 µL
Measuring range	0.05 - 200 ng/mL
Calibrators & Controls frequency	Every 28 days

IMPORTANT INFORMATION

VIDAS® B·R·A·H·M·S PCT™ (PCT) is not indicated to be used as a standalone diagnostic assay and should be used in conjunction with clinical signs and symptoms of infection and other diagnostic evidence. Decisions regarding antibiotic therapy should NOT be based solely on procalcitonin concentrations. PCT results should always be interpreted in the context of the clinical status of the patient and other laboratory results.

REFERENCES

1. Christ-Crain M, *et al.* Lancet 2004; 363:600-607.
2. Fleming-Dutra KE, *et al.* JAMA. 2016 May 3;315(17):1864-73.
3. Antibiotics – side effects. NHS Choices. <http://www.nhs.uk/Conditions/Antibiotics-penicillins/Pages/Side-effects.aspx>. Accessed February 15, 2017.
4. Lessa FC, *et al.* N Engl J Med. 2015; 372:825-834.
5. Kumar A, *et al.* Crit Care Med. 2006;34(6):1589-1596.
6. Müller B, *et al.* J Clin Endocrinol Metab. 2001;86(1):396-404.
7. Harbarth S, *et al.* Am J Respir Crit Care Med 2001;164(3):394-402.
8. Müller B, *et al.* Crit Care Med. 2000;28(4):977-983.
9. Brunkhorst, *et al.* Intensive Care Med. 1998;24(8):888-889.
10. Meisner M. J Lab Med. 1999;23:263-272.
11. Schuetz P, *et al.* Arch of Int Medicine 2011;171(15):1322-1331.
12. Albrich WC, *et al.* Arch Intern Med. 2012;172(9):715-722.
13. Schuetz P, *et al.* Cochrane Database Syst Rev. 2017;10:CD007498. doi: 10.1002/14651858.CD007498.
14. Schuetz P, *et al.* Lancet Infect Dis. 2017 S1473-3099(17)30592-3.
15. Schuetz P, *et al.* Clin Chem Lab Med. 2015 Mar;53(4):583-92.
16. Stojanovic I, *et al.* Clin Chem Lab Med. 2016 Sep 22. doi: 10.1515/cclm-2016-0349.
17. Van Der Maas *et al.* Omics: Journal of Integrated Biology Vol21, Number 4, 2017.
18. Schneider JE *et al.* Value in Health 19 (2016) A1-A318.